

*w*<sup>b</sup>

---

*b*

**UNIVERSITY  
OF BERN**

# *u<sup>b</sup>* X-ray microtomography

## 9256-HS2024-0: Advanced Microscopy

**David Haberthür**

Institute of Anatomy, December 20, 2024

# *u<sup>b</sup>* Grüessech mitenang!

- David Haberthür
  - Physicist by trade
  - PhD in high resolution imaging of the lung, Institute of Anatomy, University of Bern, Switzerland
  - Post-Doc I: TOMCAT, Swiss Light Source, Paul Scherrer Institute, Switzerland
  - Post-Doc II: µCT group, Institute of Anatomy, University of Bern, Switzerland

$u^b$


# Grüessech from the $\mu$ CT group



[David.Haberthuer@unibe.ch](mailto:David.Haberthuer@unibe.ch) [Ruslan.Hlushchuk@unibe.ch](mailto:Ruslan.Hlushchuk@unibe.ch) [Oleksiy.Khoma@unibe.ch](mailto:Oleksiy.Khoma@unibe.ch)

# $\mu$ CT-group

- microangioCT [1]
  - Angiogenesis: heart, musculature [2] and bones
  - Vasculature: (mouse) brain [3], (human) nerve scaffolds [4], (human) skin flaps [5] and tumors
- Zebrafish musculature and gills [6]
- (Lung) tumor detection and metastasis classification [7]
- Collaborations with museums [8] and scientist at UniBe [9], [10] to scan a wide range of specimens, from human hearing bones to meteorites
- Automate all the things! [11], [12]

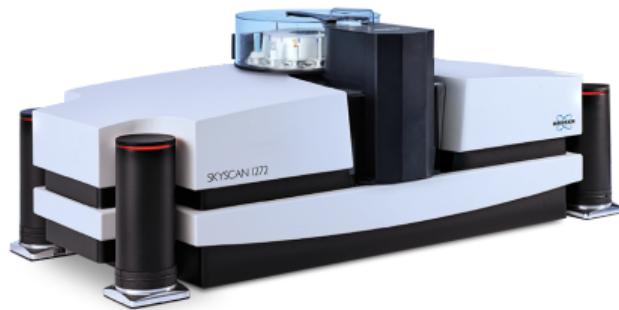

[bruker.com/skyscan1272](http://bruker.com/skyscan1272)

*u*<sup>b</sup>

# Contents

Overview

  Imaging methods

Tomography

  History

  Tomography today

A scan, from *getting started* to *nice image*

Example: A study about teeth

  Overview


  Materials & Methods

  Results

*u<sup>b</sup>*

# micro-Computed tomography

- Allows for imaging dense and non-transparent samples
- Non-destructive imaging
- Results in three-dimensional images
- Covers a very large range of sample sizes

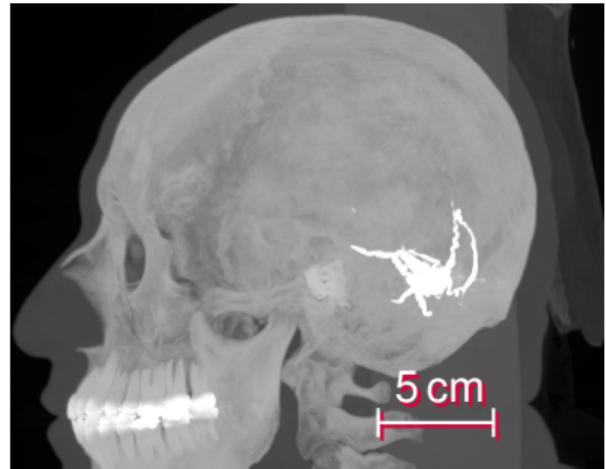


From [13], Subject C3L-02465

*u<sup>b</sup>*


# micro-Computed tomography

- Allows for imaging dense and non-transparent samples
- Non-destructive imaging, thus compatible with routine sample preparation
- Results in three-dimensional images with  $\mu\text{m}$  resolution
- Covers a very large range of sample sizes
- (Small) biological samples
- Enables correlative imaging pipelines, scanning of precious biological samples, as well as museum & collection material

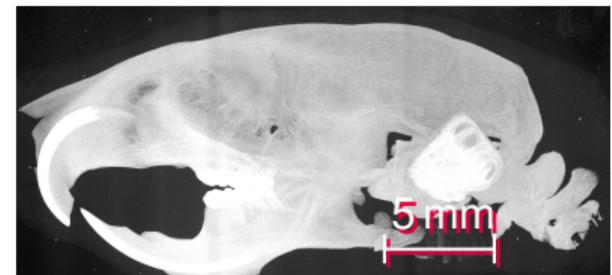

*u*<sup>b</sup>

# Imaging methods

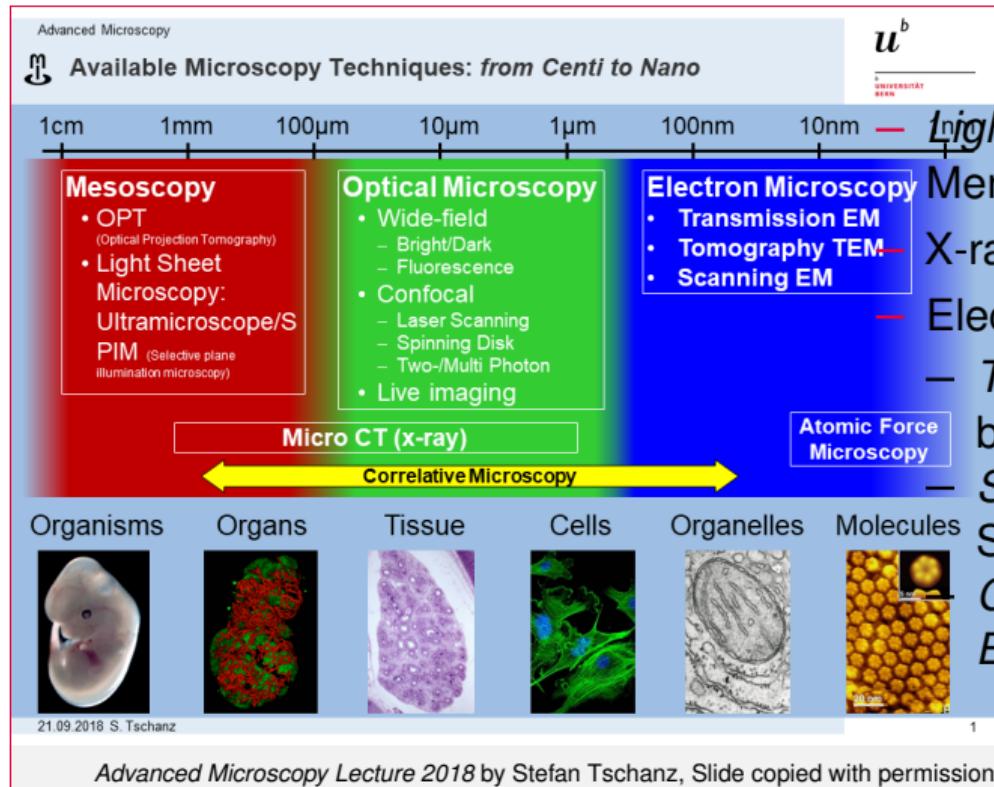

- Light Sheet Microscopy* by Nadia Mercader Huber
- X-ray imaging
- Electron microscopy
- *Transmission Electron Microscopy* by Dimitri Vanhecke
- *Scanning Electron Microscopy* by Sabine Kässmeyer and Ivana Jaric
- Cryoelectron Microscopy & Serial Block Face SEM* by Ioan Iacovache

$u^b$

# CT-Scanner



[youtu.be/2CWpZKuy-NE](https://youtu.be/2CWpZKuy-NE)

*u<sup>b</sup>*

# CT History


- 1895: Wilhelm Conrad Röntgen discovers X-rays



w.wiki/BHAN

# CT History


- 1895: Wilhelm Conrad Röntgen discovers X-rays
- **Cormack1963**: Cormack1963 used a collimated  $^{60}\text{Co}$  source and a Geiger counter as a detector **Cormack1963**
- **Hounsfield1976a**: Hounsfield1976a worked on first clinical scanner **Hounsfield1976a**

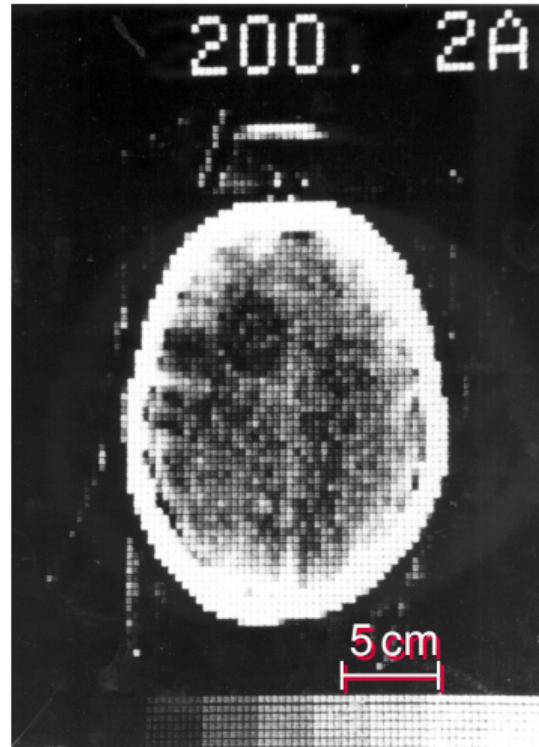

From Beckmann2006, Figure 5

# CT History


- 1895: Wilhelm Conrad Röntgen discovers X-rays
- Cormack1963: Cormack1963 used a collimated  $^{60}\text{Co}$  source and a Geiger counter as a detector Cormack1963
- Hounsfield1976a: Hounsfield1976a worked on first clinical scanner Hounsfield1976a
- CT scanner generations
  - First generation

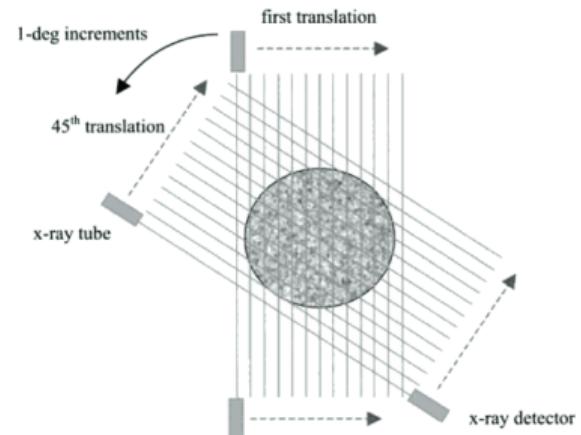

From Hsieh2003, Figure 1.12

# CT History


- 1895: Wilhelm Conrad Röntgen discovers X-rays
- **Cormack1963**: Cormack1963 used a collimated  $^{60}\text{Co}$  source and a Geiger counter as a detector **Cormack1963**
- **Hounsfield1976a**: Hounsfield1976a worked on first clinical scanner **Hounsfield1976a**
- CT scanner generations
  - First generation
  - Second generation

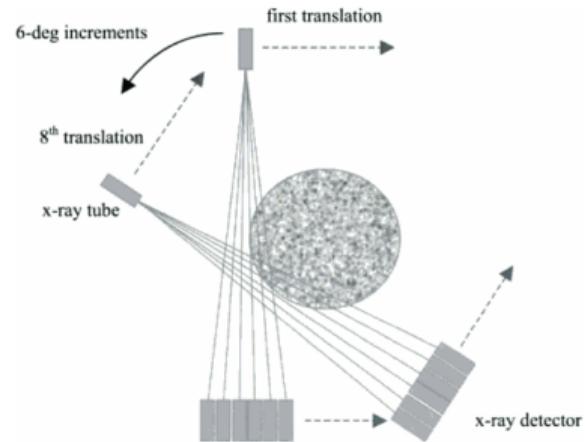

From **Hsieh2003**, Figure 1.13

# CT History

- 1895: Wilhelm Conrad Röntgen discovers X-rays
- Cormack1963: Cormack1963 used a collimated  $^{60}\text{Co}$  source and a Geiger counter as a detector Cormack1963
- Hounsfield1976a: Hounsfield1976a worked on first clinical scanner Hounsfield1976a
- CT scanner generations
  - First generation
  - Second generation
  - Third generation

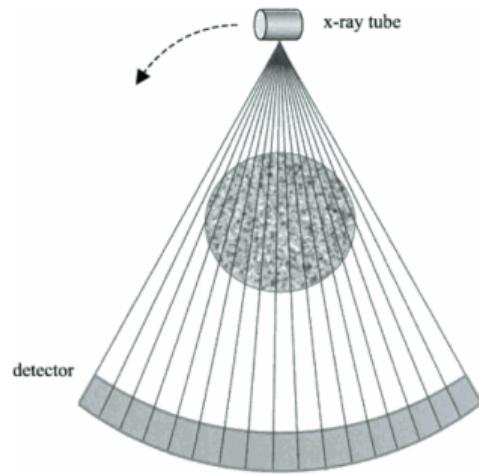

From Hsieh2003, Figure 1.14

# $\mu$ CT History I


- X-ray computed tomography began to replace analog focal plane tomography in the early 1970s **Lin2019**
- Non-medical use in the late 1970s, for detection of internal defects in fabricated parts and equipment
- Lee Feldkamp **Feldkamp1984** developed one of the early laboratory microCT systems by assembling a micro-focus cone beam x-ray source, specimen holder and stages, and an image intensifier at Ford Motor Company's Scientific Research Laboratory to nondestructively detect damage in ceramic manufactured automobile parts
- Feldkamp met with scientists at Henry Ford Hospital and University of Michigan interested in understanding the relationship between the microstructure and biomechanical function of trabecular bone to study osteoporotic fractures **Feldkamp1983**

# $\mu$ CT History II


- $\mu$ CT was first reported in the 1980s, for scanning gemstones to cut out the largest possible one
- Today: Nondestructive imaging for quantifying the (micro)structure of organic materials
  - Mineralized bone tissue and the relationships between the mechanical behavior of bone to its structural and compositional properties
  - Teeth and their internal details
  - Soft tissues and vasculature using radio-opaque contrast agents
  - Characterization of anatomical details in high resolution
- $\approx$ 2500  $\mu$ CT systems are in use worldwide with over 1000 publications annually

$u^b$

# Why $\mu$ CT?

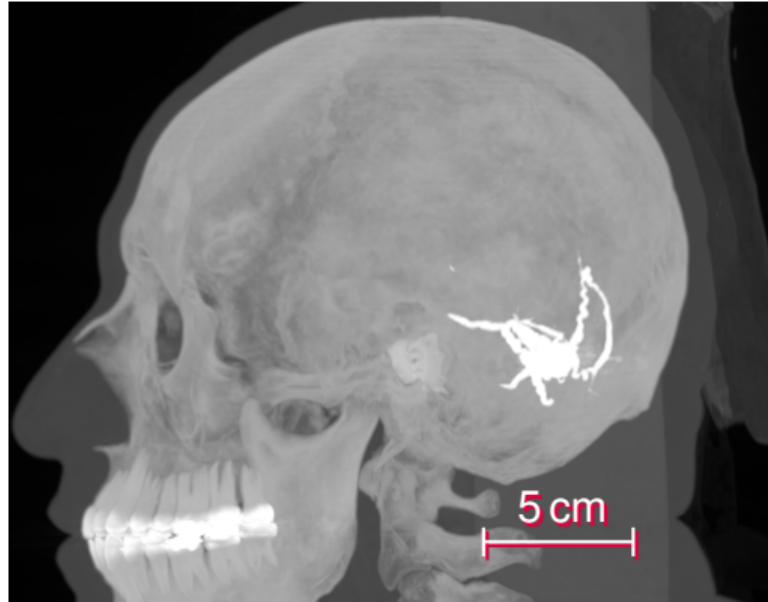


From [13], Subject C3L-02465

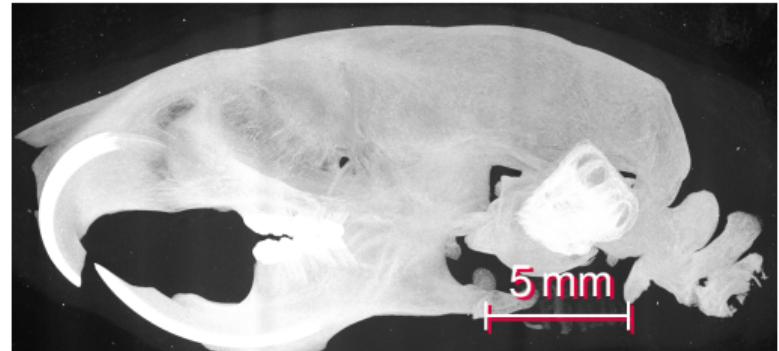

$u^b$

# Why $\mu$ CT?

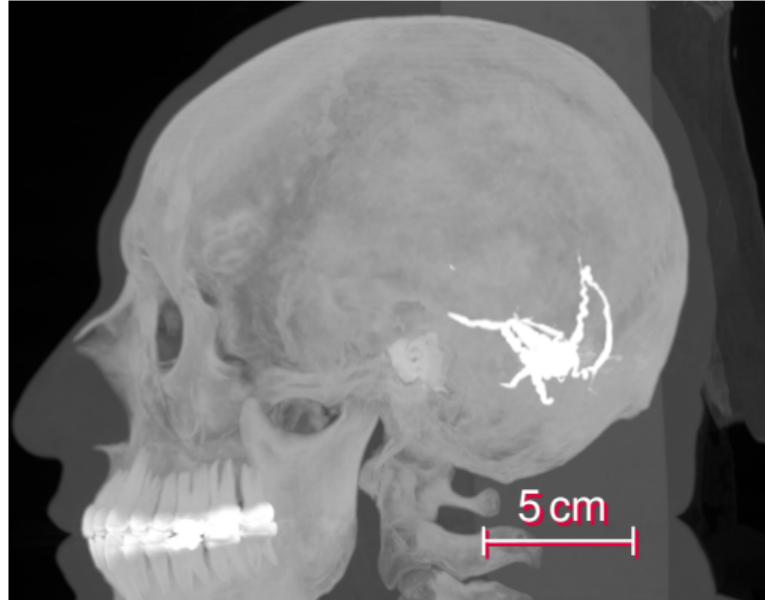


From [13], Subject C3L-02465




$u^b$

# Why $\mu$ CT?

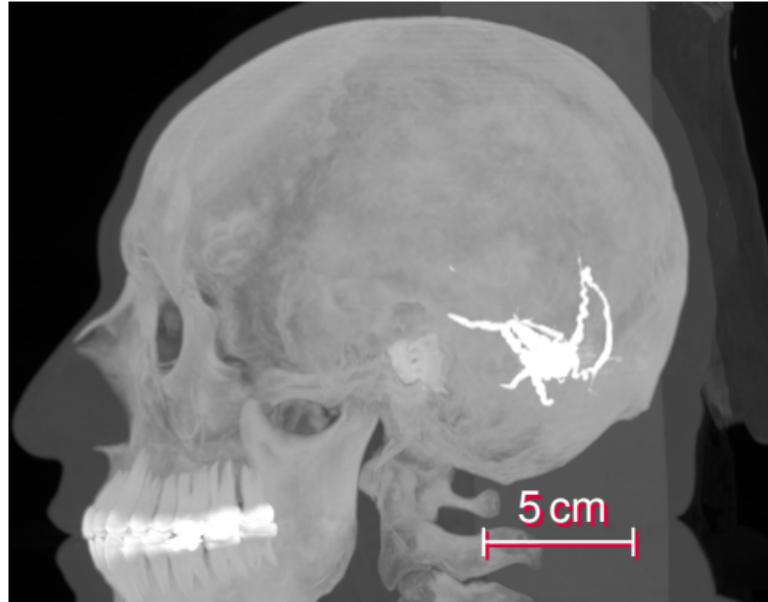


From [13], Subject C3L-02465

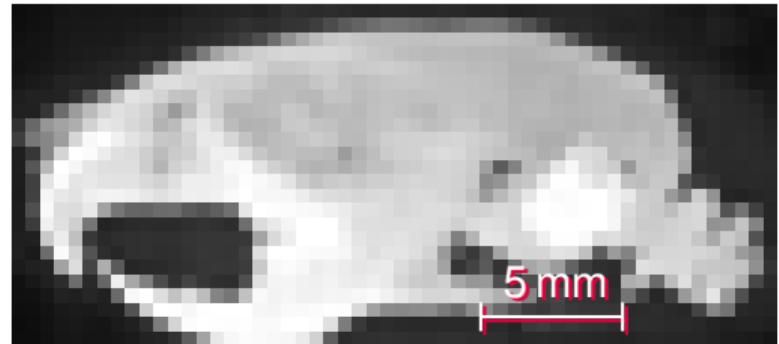


$u^b$

# Why $\mu$ CT?

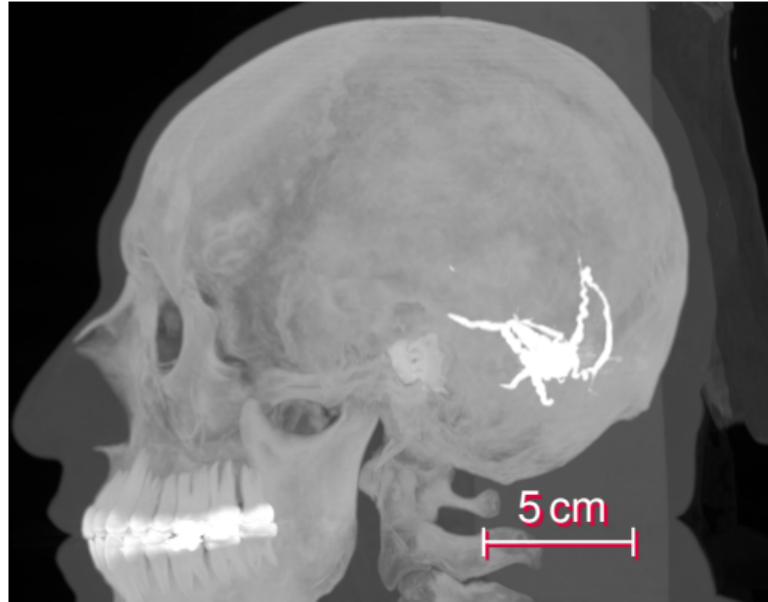


From [13], Subject C3L-02465

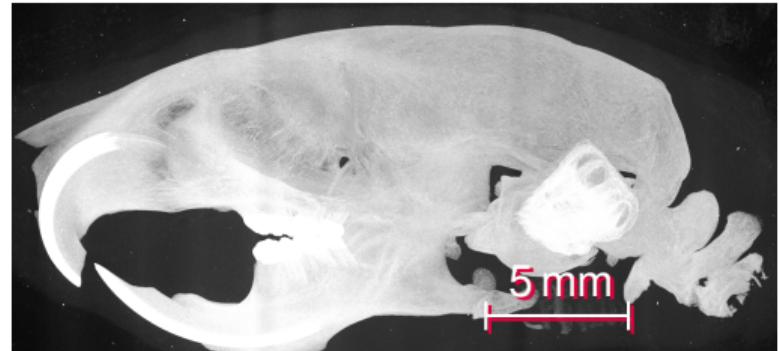


$u^b$

# Why $\mu$ CT?

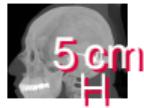




From [13], Subject C3L-02465

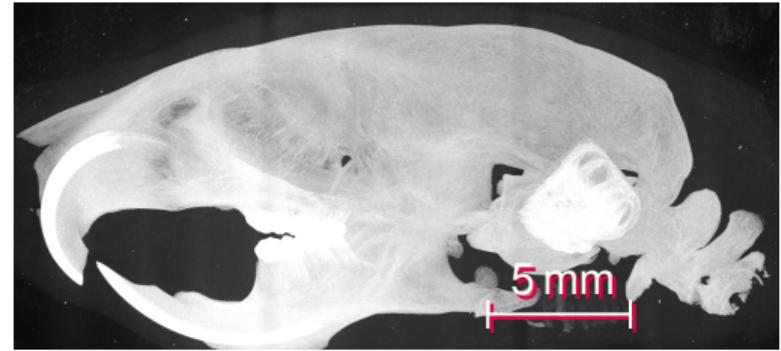

# Maximum intensity projection

## Projections

Reducing the dimensions of a dataset. For example projecting a volume (3D) to a surface by taking the maximum value across planes for each pixel.

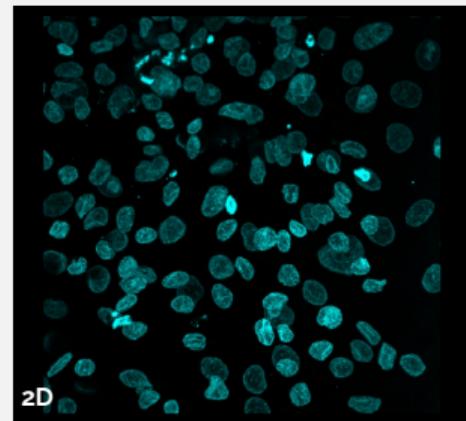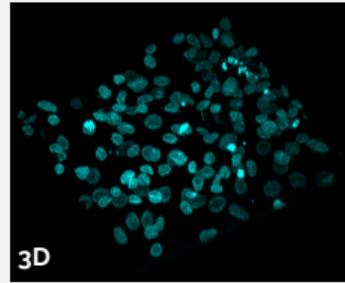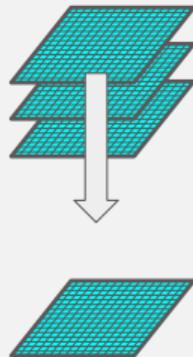

# Machinery

- Hospital CT
  - Voxel size around 0.5 mm
- Lab/Desktop CT
  - Voxel size around 7  $\mu\text{m}$  (*in vivo*)
  - Voxel size around 0.5  $\mu\text{m}$  (*ex vivo*)
- Synchrotron CT
  - Voxel size down to 160 nm



flic.kr/p/D4rbom

# Machinery


- Hospital CT
  - Voxel size around 0.5 mm
- Lab/Desktop CT
  - Voxel size around 7  $\mu\text{m}$  (*in vivo*)
  - Voxel size around 0.5  $\mu\text{m}$  (*ex vivo*)
- Synchrotron CT
  - Voxel size down to 160 nm



[flic.kr/p/fpTrGu](https://flic.kr/p/fpTrGu)

# Machinery

- Hospital CT
  - Voxel size around 0.5 mm
- Lab/Desktop CT
  - Voxel size around 7  $\mu\text{m}$  (*in vivo*)
  - Voxel size around 0.5  $\mu\text{m}$  (*ex vivo*)
- Synchrotron CT
  - Voxel size down to 160 nm

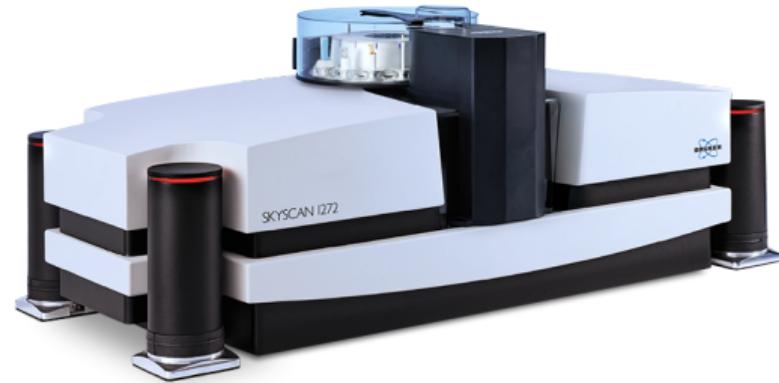

[bruker.com/skyscan1272](http://bruker.com/skyscan1272)

# Machinery

- Hospital CT
  - Voxel size around 0.5 mm
- Lab/Desktop CT
  - Voxel size around 7  $\mu\text{m}$  (*in vivo*)
  - Voxel size around 0.5  $\mu\text{m}$  (*ex vivo*)
- Synchrotron CT
  - Voxel size down to 160 nm




flic.kr/p/7Xhk2Y © ⓘ ⓘ

$u^b$

# What is happening?

No matter what kind of machine, the basic principle is always

- an x-ray source
- a sample
- a detector

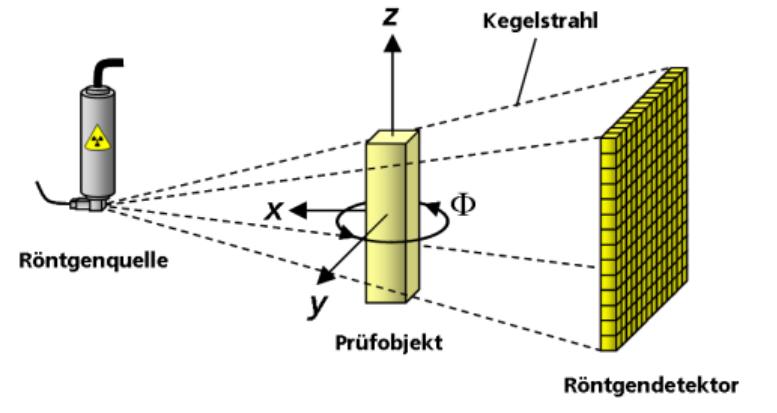


# Magnification

## Introduction – why do we need microscopes?

4

- Basic geometric optics: how does the eye form an image?

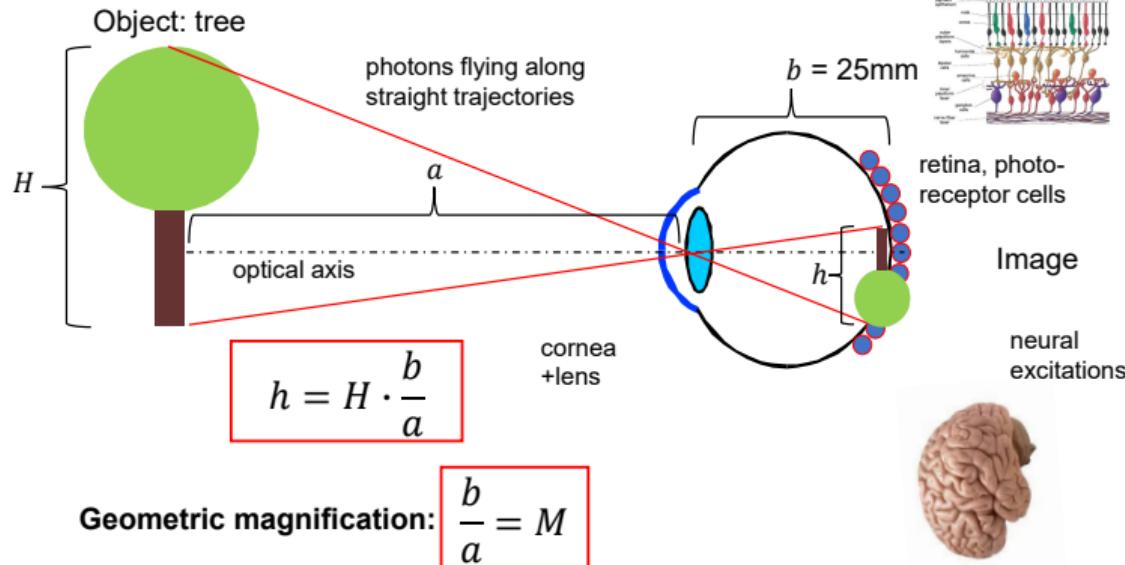

*u<sup>b</sup>*


# Preparation

- Study design
- Sample preparation

# Projections


## X-ray generation and contrast

X-ray tube:  
nearly point like  
photon source

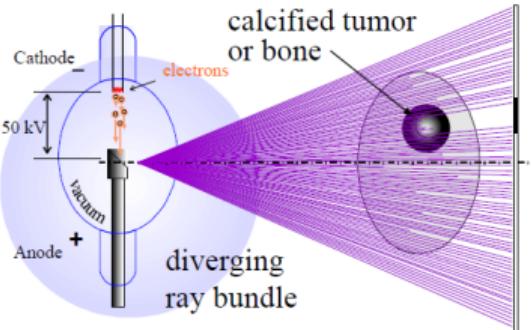


Contrast is given by  
absorption of intensity I

Note that contrast is negative  
X-ray shadowgraphy  
is a bright field technique

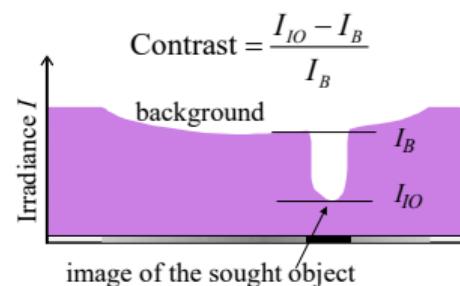


# $u^b$ Projection acquisition

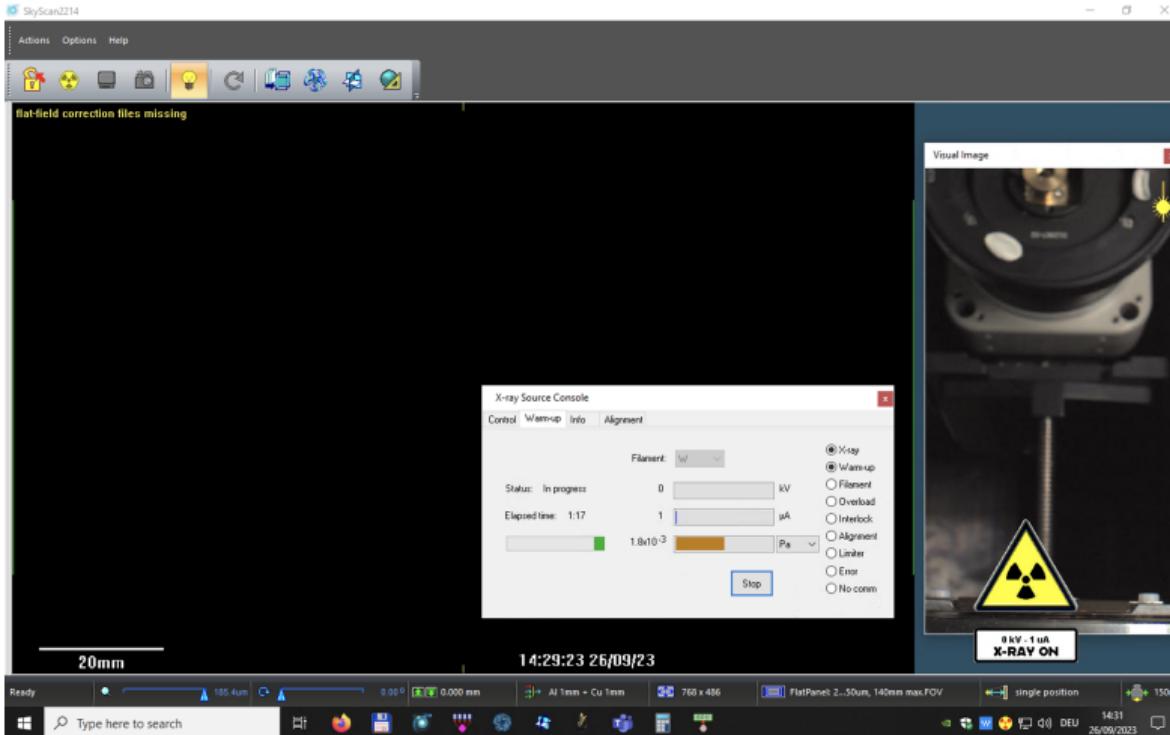

$u^b$

# Projections

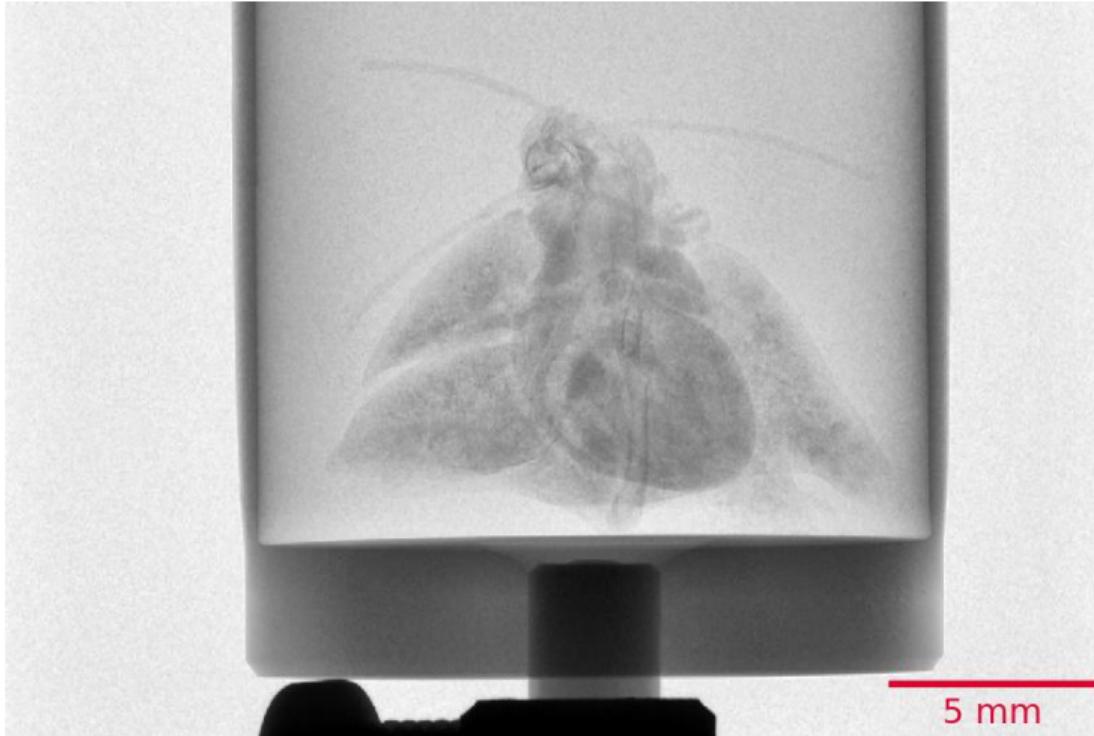


# Projections

- A (micro-focus) x-ray source illuminates the object
- The x-rays penetrate the sample and are attenuated
- A scintillator converts the x-rays to visible light
- A (planar) x-ray detector collects (magnified) projection images.
- The projections are recorded on disk


*u*<sup>b</sup>

# Reconstructions

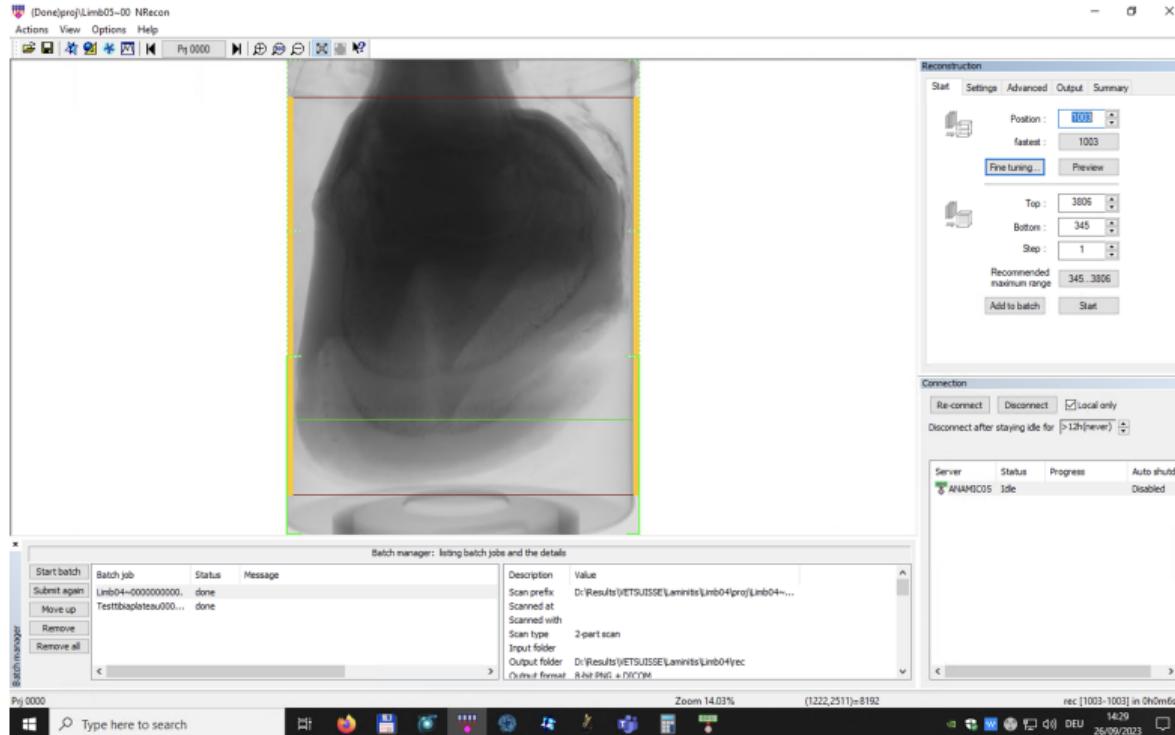


*u*<sup>b</sup>

# Reconstructions

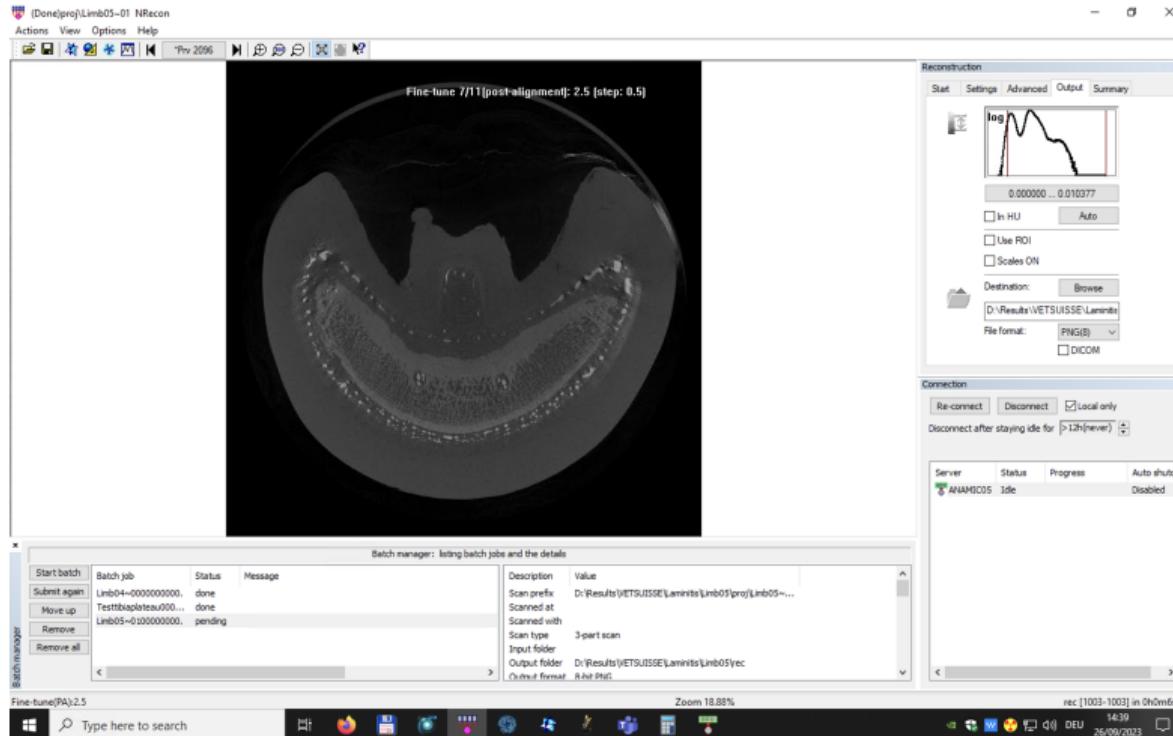

*u*<sup>b</sup>

# Reconstructions

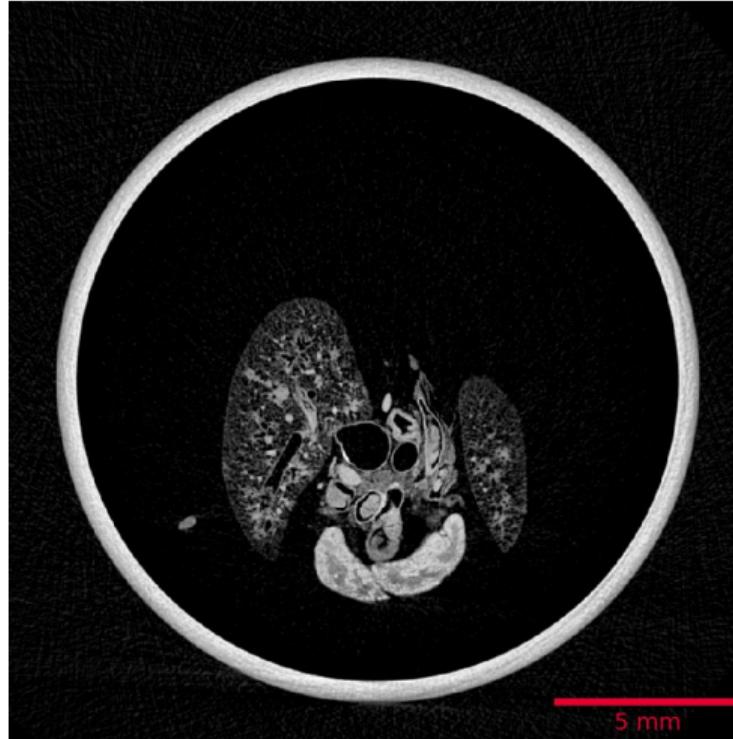

# Reconstructions

- Based on hundreds of angular views acquired while the object rotates, a computer synthesizes a stack of virtual cross section slices through the object.
- Radon Transformation
- Filtered back projection
- Cone beam reconstruction **Feldkamp1984**
- Corrections (beam hardening, etc.)
- Writing to stack

$u^b$

# Visualization



# Visualization


- Based on reconstructions, a computer synthesizes a three-dimensional view of the scanned sample

$u^b$

# Imaging

CEM 2024-Light Sheet Microscopy

## How to deal with large multidimensional data?

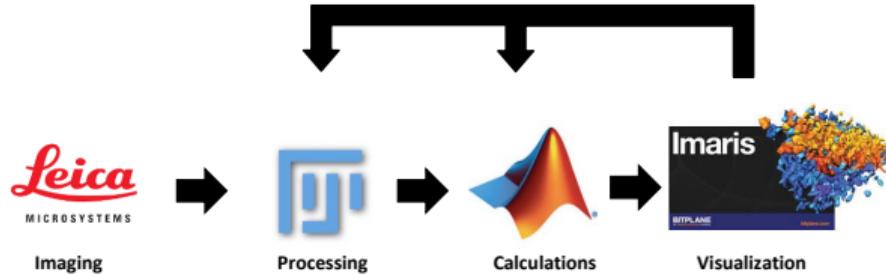

Customized macro Beatsync (Michael Liebling)  
for file conversion Heart synchronization

CEM 2024-Light Sheet Microscopy

*Light Sheet Microscopy by Nadia Mercader, Slide 43*

# What to use?

- ImageJ/Fiji **Schindelin2012**
- Also see *Fundamentals of Digital Image Processing* by Guillaume Witz
- Reproducible research
  -  in Jupyter **Kluyver2016**
  - **git**
  - Script all your things!
  - Data repositories; i. e. sharing is caring!

# Quantitative data

- Pretty images are nice, but we need quantitative numbers
- Segmentation
- Characterization

*u<sup>b</sup>*

# Internal morphology of human teeth

Collaboration with zmk bern – Zahnmedizinische Kliniken

- Number of teeth
- Structure of teeth



- Morphology of teeth
- Analytical workflow, [doi.org/g/r8](https://doi.org/g/r8)


# How?

- 104 extracted human permanent mandibular canines
- $\mu$ CT imaging
- Root canal configuration, according to **Briseno-Marroquin2015** Briseno-I
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!



# How?

- 104 extracted human permanent mandibular canines
- $\mu$ CT imaging
- Root canal configuration, according to **Briseno-Marroquin2015 Briseno-M**
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!

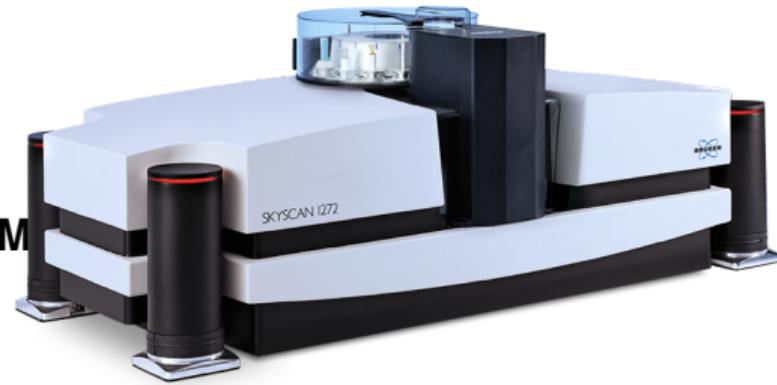

[bruker.com/skyscan1272](http://bruker.com/skyscan1272)

- 104 extracted human permanent mandibular canines
- µCT imaging
- Root canal configuration, according to **Briseno-Marroquin2015** **Briseno-Marroquin2015**
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!


```
Scanner=SkyScan1272
Instrument S/N=15G09089-B
Software Version=1.1.19
Filename Prefix=Tooth045~00
Number Of Files= 482
Number Of Rows= 1092
Number Of Columns= 1632
Source Voltage (kV)= 80
Rotation Speed (deg/s)= 125
Image Pixel Size (µm)=9.999986
Exposure (ms)=950
Rotation Step (deg)=0.400
Frame Averaging=ON (3)
Filter=Al 1mm
Study Date and Time=02 Jul 2020 08h:23m:34s
Scan duration=0h:39m:51s
```

# How?

- 104 extracted human permanent mandibular canines
  - $\mu$ CT imaging
  - Root canal configuration, according to **Briseno-Marroquin2015** ~~Briseno-Marroquin2015~~
  - *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!
- Sample changer on the SkyScan 1272*  
In total:
- 13 days of *continuous*  $\mu$ CT scanning
  - 819 GB of raw data
  - 230 648 TIFF projections
  - 326 GB data as input for analysis
  - 282 062 PNG reconstructions

# How?


- 104 extracted human permanent mandibular canines
- µCT imaging
- Root canal configuration, according to **Briseno-Marroquin2015 Briseno-I**
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!

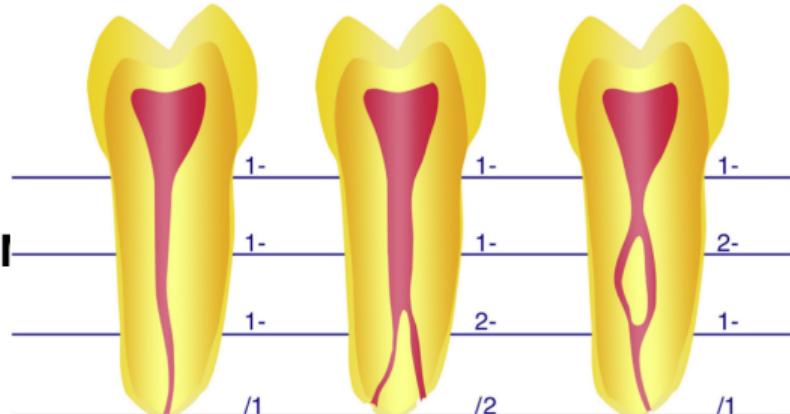

From **Briseno-Marroquin2015**, Fig. 2

# How?

- 104 extracted human permanent mandibular canines
- $\mu$ CT imaging
- Root canal configuration, according to **Briseno-Marroquin2015 Briseno-I**
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!

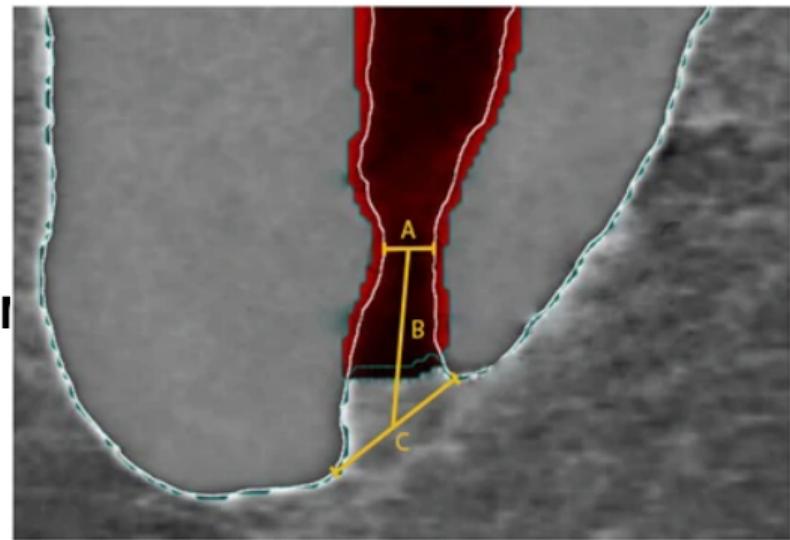

From **Wolf2017**, Fig. 1

*u*<sup>b</sup>

# How?

- 104 extracted human permanent mandibular canines
- µCT imaging
- Root canal configuration, according to **Briseno-Marroquin2015** Briseno-I
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!



[gph.is/2nqkple](https://gph.is/2nqkple)

*u*<sup>b</sup>

# How?

- 104 extracted human permanent mandibular canines
- $\mu$ CT imaging
- Root canal configuration, according to **Briseno-Marroquin2015** Briseno-M
- *Reproducible* analysis **Haberthuer2020a**, e.g. you can click a button to double-check or recalculate the results yourself!

The screenshot shows a GitHub repository interface. At the top, there are buttons for 'master', '1 branch', '1 tag', 'Go to file', 'Add file', and a green 'Code' button. Below this is a list of commits from a user named 'habi'. The commits are as follows:

File	Commit Message	Time Ago
.github/actionsfile	Update actions file	20 days ago
.gitignore	Only 'mode' changes	2 months ago
DownloadFromOSF.ipynb	Clean run of download script	22 days ago
README.md	Typo in Binder badge & link to full repo on Binder	22 days ago
Tooth.Border.jpg	Only 'mode' changes	2 months ago
Tooth.Characterization.jpg	Only 'mode' changes	2 months ago
ToothAnalysis.ipynb	Only select a subset if we actually have data:wq	22 days ago
ToothDataSize.ipynb	Clean run of notebook	22 days ago
ToothDisplay.ipynb	Display Tooth045 for manuscript	22 days ago
requirements.txt	We also need this	2 months ago
treebeard.yaml	Add treebeard configuration	20 days ago


Below the commits is a 'README.md' file. It contains a DOI link (10.5281/zenodo.3999402), a 'treebeard.yml' status indicator (failing), and a 'Launch binder' button. The file also includes a section titled 'A big tooth cohort' with the following text:

We scanned a big bunch of teeth for a team of the dental clinic of the University of Bern.

To get an overview of the samples while we scanned the whole tooth cohort we generated a [preview](#) and [analysis notebook](#). The analysis notebook (with download possibility for two of the +100 teeth) can be started in your browser by clicking on the 'Binder' badge above, without installation of any software. If you'd like to start a Binder instance with the full repository, you can click [here](#).

$\mu$ b



# $\mu$ CT imaging

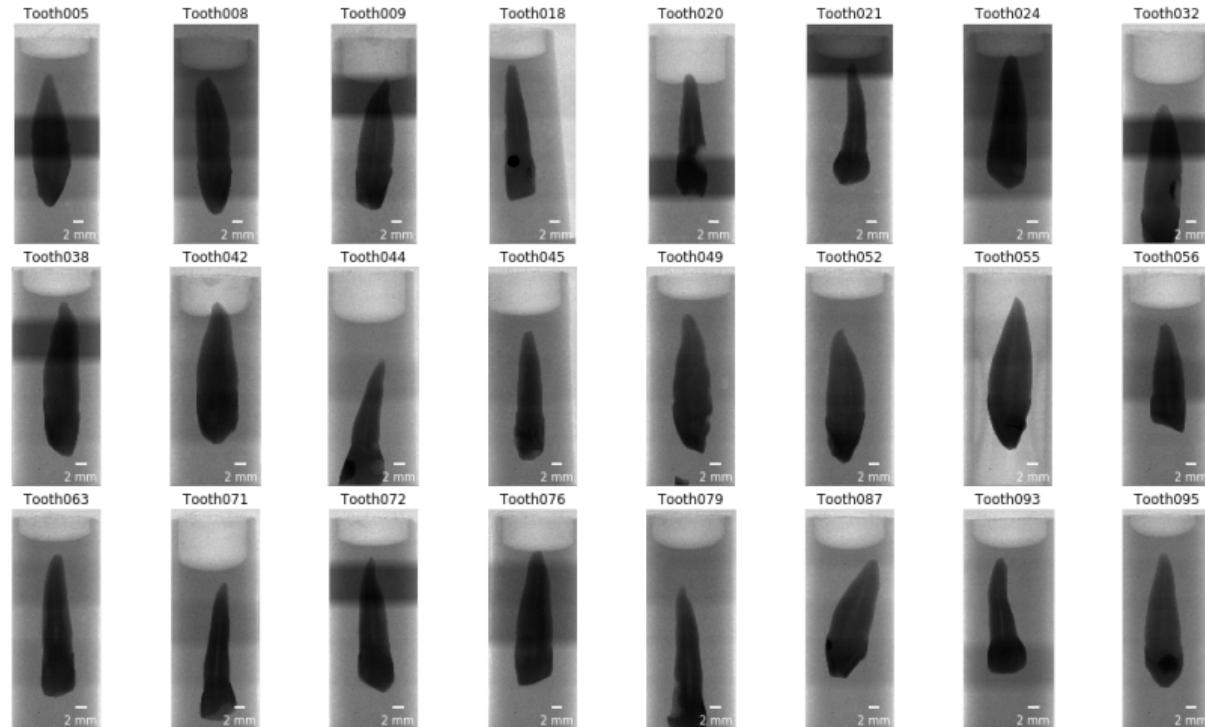

*u*<sup>b</sup>

# Dataset cropping




- Full datasets: 326 GB
- Cropped datasets: 115 GB

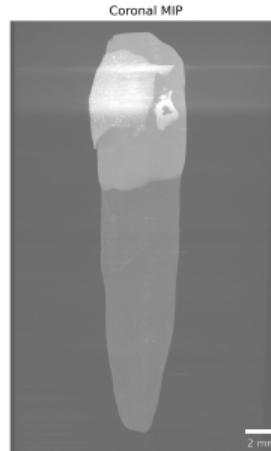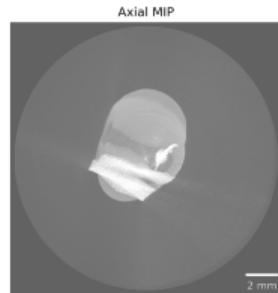

*u*<sup>b</sup>

# Dataset cropping

- Full datasets: 326 GB
- Cropped datasets: 115 GB

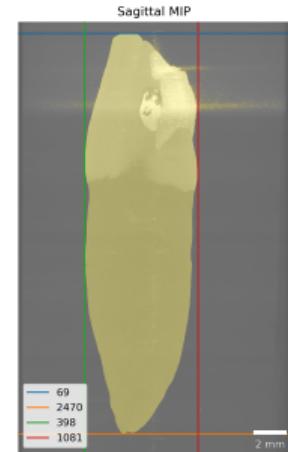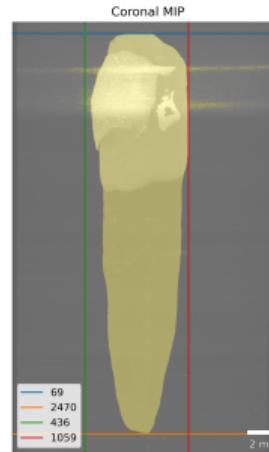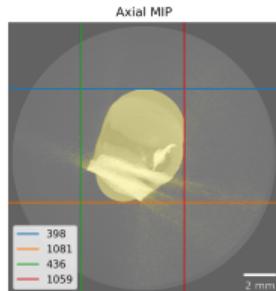


*u*<sup>b</sup>

# Tooth morphology



*u<sup>b</sup>*


# Tooth morphology

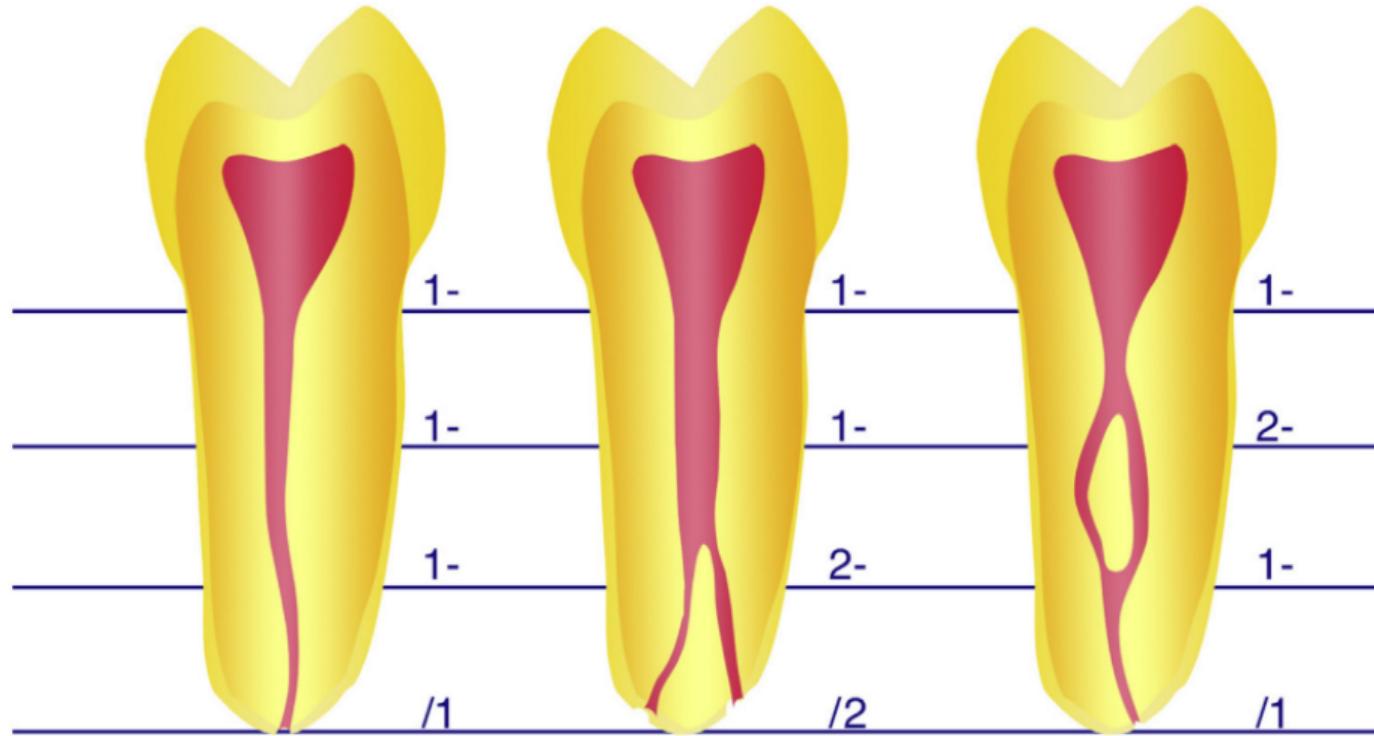

From Briseno-Marroquin 2015, Fig. 2




*u<sup>b</sup>*

# Tooth morphology

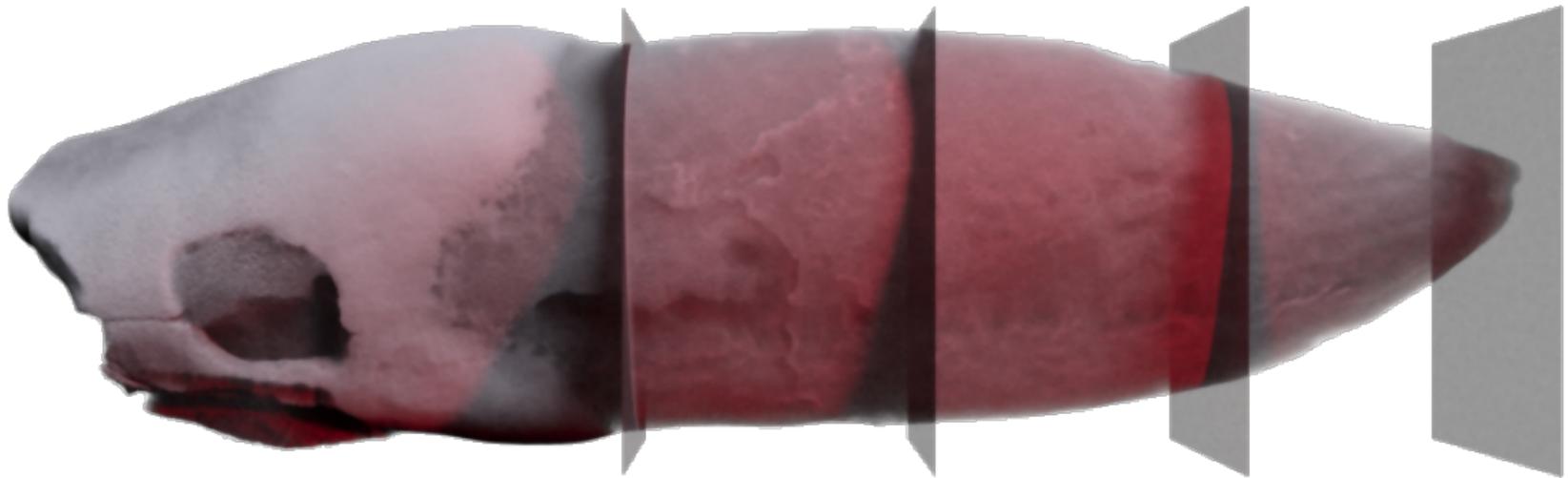


*u<sup>b</sup>*

# Detection of enamel-dentin border

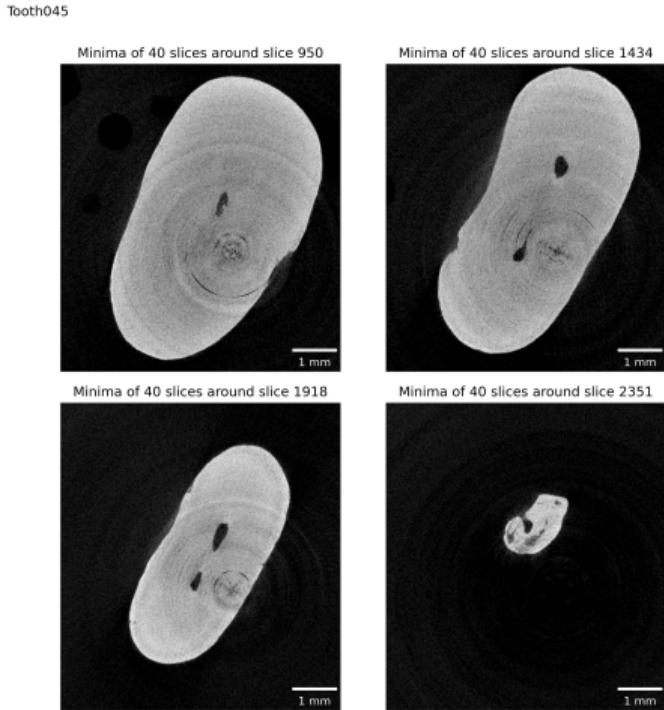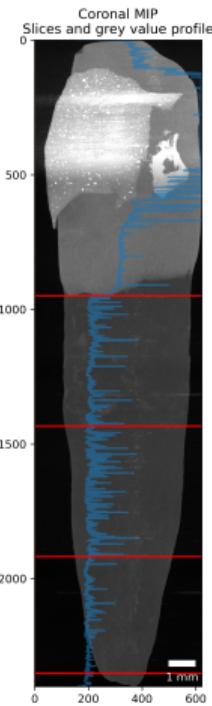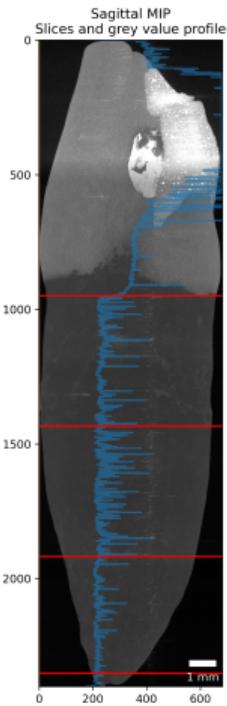

*u*<sup>b</sup>

# Detection of enamel-dentin border

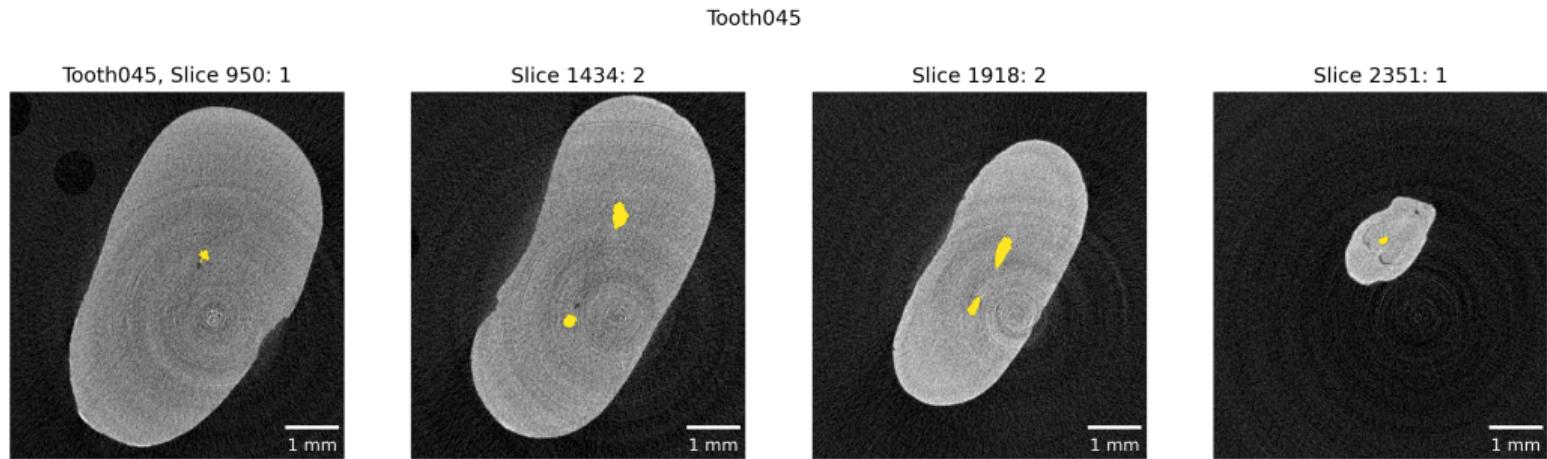


*u<sup>b</sup>*

# Outcome root canal configuration classification

Roots	RCC	#	%	
Single (N=98)	1-1-1/1	73	74.5	
	1-1-1/2	14	14.3	
	1-1-1/3	1	1.0	
	1-1-1/4	2	2.1	
	1-1-2/1	1	1.0	
	1-2-1/1	4	4.1	
	1-2-1/2	1	1.0	
	1-2-2/2	1	1.0	
	2-3-1/1	1	1.0	
Double (N=3)	Buccal	1-1-1/1	2	66.6
		1-2-1/1	1	33.3
	Lingual	1-1-1/1	2	66.6
		1-1-1/2	1	33.3

*u*<sup>b</sup>

# Extraction of root canal space

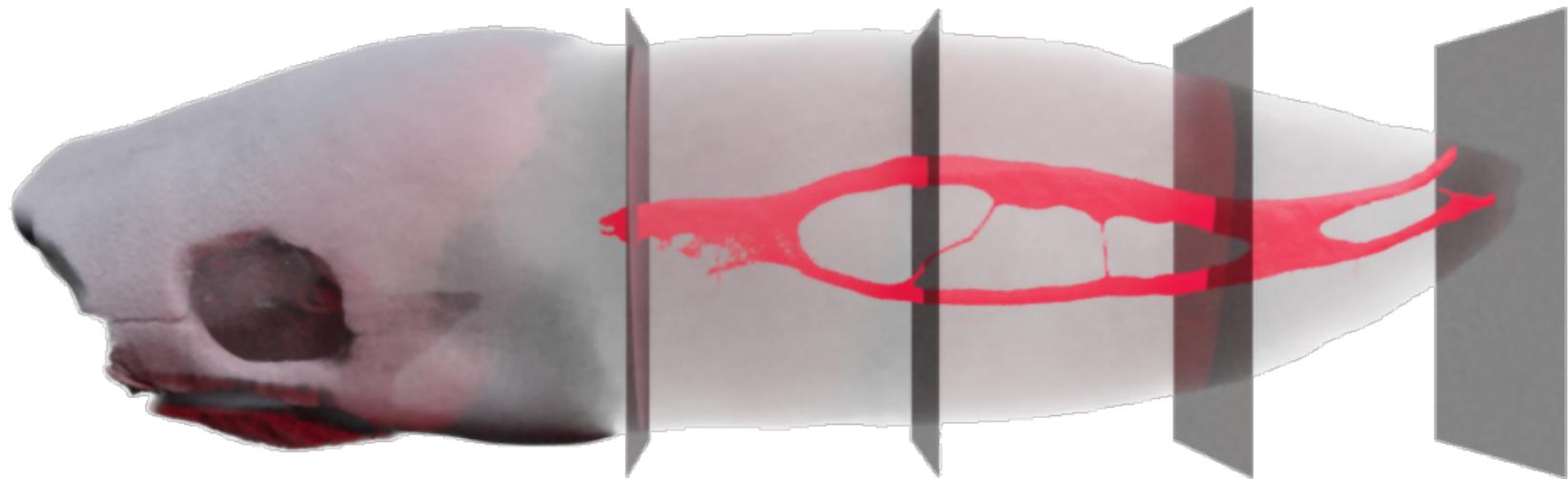

# Conclusion ZMK

- Efficient use of time, e. g. more teeth does not mean more (human) work
- Reproducible analysis with *free and open-source* software, usable by *anyone*
- Objective analysis, e. g. no operator bias

*u*<sup>b</sup>

# Thanks!

- Thanks for listening to me!
- What questions do you have for me?

# Colophon

- This BEAMER presentation was crafted in  $\text{\LaTeX}$  with the (slightly adapted) template from *Corporate Design und Vorlagen* of the University of Bern.
  - Complete source code: [git.io/fjpP7](https://git.io/fjpP7)
  - The  $\text{\LaTeX}$  code is automatically compiled with a GitHub action to a (handout) PDF which you can access here: [git.io/JeQxO](https://git.io/JeQxO)
- Did you spot an error?
  - File an issue: [git.io/fjpPb](https://git.io/fjpPb)
  - Submit a pull request: [git.io/fjpPN](https://git.io/fjpPN)
  - Send me an email: [david.haberthuer@unibe.ch](mailto:david.haberthuer@unibe.ch)